JOHN BUNN®  
A Graham-Field Brand  
DigiOx  
Finger Pulse Oximeter

Operation Manual

**Important: Do not operate the John Bunn DigiOx Finger Pulse Oximeter without first reading and understanding this manual! Save this manual for future use.**

**Note: The most current version of this manual can be found online at [www.grahamfield.com](http://www.grahamfield.com)**

JB02008-INS-LAB-RevA09

▲ Cautions

To ensure patient safety, use only parts and accessories specified in this manual.

Keep the oximeter away from dust, vibration, corrosive substances, explosive materials, high temperature and moisture.

If the oximeter gets wet, stop operation immediately and do not resume operation until it is dry. When it is carried from a cold environment to a warm and/or humid environment, allow it to warm to room temperature before using it. DO NOT use sharp objects to press the front panel button.

DO NOT sterilize this device. Do not subject this device to irradiation, steam, ethylene oxide, high temperature, or immersion in liquid. Refer to the Repair and Maintenance/Cleaning and Disinfection section of this manual for instructions for cleaning and disinfection.

6. To change the oximeter's display screen brightness, press the Power/Display Button for more than one second. There are 10 brightness levels, and the level, from "BR 1" to "BR 10", will show on the screen during brightness adjustment. The default is level four. Keep the Power/Display Button depressed for as long as it takes to cycle to the desired brightness level.

7. The oximeter will automatically turn off if no finger is sensed for 8 seconds. If the finger is removed and the Power/Display Button is pressed, the oximeter will remain on and change screens, but will display "FINGER OUT".

ACCESSORIES (INCLUDED)

- 1 One lanyard
- 2 Two AAA alkaline batteries
- 3 One operation manual

REPAIR, MAINTENANCE, AND DISPOSAL

⚠ **WARNING:** DO NOT open the oximeter housing. It contains no serviceable parts. Service must be performed by trained, authorized personnel only. Otherwise, device failure and health hazard may occur.

Maintenance


In order to ensure the oximeter's longevity and continued dependable operation, please observe the following maintenance recommendations.

CONTENTS

INTRODUCTION .....	3
INTENDED USE .....	3
CONTRAINDICATIONS.....	3
IMPORTANT SAFETY GUIDELINES - PLEASE READ BEFORE USE .....	3
OVERVIEW .....	7
APPEARANCE .....	7
FEATURES .....	7
NAME AND MODEL .....	8
BATTERY INSTALLATION .....	8
OPERATION INSTRUCTIONS .....	8
ACCESSORIES (INCLUDED) .....	11
REPAIR, MAINTENANCE, AND DISPOSAL.....	11
MAINTENANCE.....	11
CLEANING AND DISINFECTION .....	12
DISPOSAL .....	13
TROUBLESHOOTING .....	13
SYMBOL KEY .....	15
TECHNICAL SPECIFICATIONS.....	16
LIMITED WARRANTY .....	18
MANUAL DE OPERACIÓN (ESPAÑOL) .....	21

OVERVIEW

Appearance

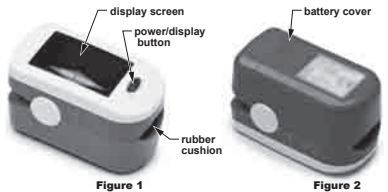The DigiOx Finger Pulse Oximeter is shown above.

Features

The John Bunn DigiOx Finger Pulse Oximeter provides a simple way to spot-check patients by combining the sensor and monitor into one integrated, compact, easy-to-use device with single-button operation. The oximeter can accurately measure oxygen saturation and pulse rate. When a finger is inserted into the sensor's rubber cushion, and the Power/Display Button is pressed to turn the oximeter on, then the SpO<sub>2</sub> value and pulse rate value display on the screen. The oximeter features a bright-color LED screen with adjustable brightness and multi-directional displays that enable the user to read results

INTRODUCTIONS

Dear User: Please read and understand all instructions before using the John Bunn DigiOx Finger Pulse Oximeter, JB02008.

Intended use

The John Bunn DigiOx Finger Pulse Oximeter is used to measure oxygen saturation (SpO<sub>2</sub>) and pulse rate in home or clinical settings. It is designed for adult and pediatric patients.

Contraindications

⚠ **WARNING:** The oximeter is not intended for continuous patient monitoring. DO NOT use the oximeter while patient is undergoing MRI or CT scanning.

Important Safety Guidelines - Please Read Before Use

The safety statements presented in this chapter refer to the basic safety information that the operator of the John Bunn DigiOx Finger Pulse Oximeter shall pay attention to and abide by. There are additional safety statements in other chapters or sections, which may be the same as or similar to the following, or specific to the operations.

⚠ **WARNING:** Indicates a potential hazard situation or unsafe practice that, if not avoided, could result in death or serious injury.

from any angle. The oximeter features a low voltage indicator, powers off automatically in eight seconds when not in use, and provides approximately 30 hours of use with the use of alkaline batteries.

Name and Model

Name: John Bunn DigiOx Finger Pulse Oximeter

Model: JB02008

BATTERY INSTALLATION

1. Refer to Figure 3; insert two AAA batteries, oriented as shown, into the battery compartment.
2. Replace the cover.

*Note: Ensure that the batteries are correctly installed. Incorrect installation will prevent the device from operating.*

OPERATION INSTRUCTIONS

Notes:

*Clean the test finger with isopropyl alcohol before and after each measurement. Ensure that the finger is dry before being inserted into the oximeter.*

*Do not shake the finger; keep finger motionless during oximeter use.*

Disposal

Do not treat this device or its batteries as household waste. Separate collection is required.

TROUBLESHOOTING

Problem	Possible Reason	Solution
The device does not turn on	SpO <sub>2</sub> and Pulse Rate display unstable or hard to see	Reposition finger correctly and try again
	Finger is shaking or patient is moving	Ensure patient remains immobile
	Patient's SpO <sub>2</sub> value is too low to measure	⚠ Remeasure; if unsuccessful, yet device is functioning correctly, contact a medical professional immediately
	Excessive ambient light	Reduce ambient light

LIMITED WARRANTY

GF Health Products, Inc. warrants the John Bunn DigiOx Finger Pulse Oximeter, Model JB02008, for a period of one year against manufacturer's defects. If a product is deemed to be under warranty, GF Health Products, Inc. shall provide, at its option, (1) replacement of any defective part or product or (2) a credit of the original selling price made to GF Health Products, Inc.'s initial customer. The warranty does not include any labor charges incurred in replacement part(s) installation or any associated freight or shipping charges to GF Health Products, Inc.

▲ **Caution:** Indicates a potential hazard or unsafe practice that, if not avoided, could result in minor personal injury or product/property damage.

*Note: Provides application recommendations or other useful information to ensure that you get the most from your product.*

⚠ WARNINGS:

**Important! Read and understand this manual before using the oximeter. Do not use this device without proper instruction from a Healthcare Professional. If the oximeter is not properly operated, personal injury and/or damage to the oximeter could result.**

If components are damaged or missing, contact your dealer immediately. DO NOT use substitute parts.

Inspect the device before use to ensure that there is no visible damage that may affect patient or user's safety or measurement performance. Thereafter, inspect the device at least once a week.

GF Health Products, Inc. assumes no responsibility for any damage or injury caused by improper installation or use of this product.

DO NOT open the oximeter housing. It contains no serviceable parts. Service must be performed by trained, authorized personnel only. Otherwise, device failure and health hazard may occur.

*Avoid placing the device on a limb which is wrapped with a cuff for blood pressure measurement or undergoing venous infusion.*

*Ensure that no obstruction (contamination, scar tissue, fingernail polish, acrylic fingernail, etc.) exists at the site where the sensor is placed; such obstruction could affect the signal received by the sensor and result in an incorrect measurement.*

*Vigorous exercise and electrosurgical device interference may affect the measuring accuracy.*

1. Open the clip as shown in Figure 4.
2. Insert finger, nail up, into the clip's rubber cushion (ensure that the finger is correctly positioned), and then gently clip the finger.
3. Press the Power/Display Button to turn the oximeter on.
4. The screen will then display the data as shown in Figure 5. After the values become stable, the user can read the information on the display screen.
5. When the screen displays as shown in Figure 5, pressing the Power/Display Button will change the screen's orientation:

**EXPLOSION HAZARD:** Do not use the oximeter in the presence of flammable anesthetics, explosive substances, vapors or liquids.

The oximeter is designed for real-time and rapid measurement. It is not suitable for long-time continuous patient monitoring. DO NOT apply the sensor to the same finger for more than two hours. If any readings are found to be abnormal, please change the position of the sensor.

The oximeter is intended only as an adjunct in patient assessment. It must be used in conjunction with clinical signs and symptoms. It is intended as a measuring, not a treatment, device.

Magnetic and electrical fields are capable of interfering with the proper performance of the oximeter. For this reason, ensure that all external devices operated in the vicinity of the oximeter comply with the relevant EMC requirements. Mobile phones, X-ray equipment, and/or MRI devices are a possible source of interference, as they may emit high levels of electromagnetic radiation. This device complies with IEC60601-1-2.

DO NOT clip this device on edema or tender tissue.

The infrared light emitted from the device (it is invisible) is harmful to the eyes. DO NOT stare into the light.

Keep this device out of the reach of children.

- a. Press the Power/Display Button once more; the screen will rotate 180° and display as shown in Figure 6.
- b. Press the Power/Display Button once more; the screen will rotate 90° and display as shown in Figure 7.
- c. Press the Power/Display Button once more; the screen will rotate 180° and display as shown in Figure 8.
- d. Press the Power/Display Button once more; the screen will rotate -90° and also show pulse as a filled wave-form bar graph along bottom of screen, as shown in Figure 9.
- e. Press the Power/Display Button once more; the screen will now show pulse as a hollow wave-form bar graph along bottom of screen, as shown in Figure 10.
- f. Press the Power/Display Button once more; the screen will return to the original screen as shown in Figure 5 on previous page.

SYMBOL KEY

Symbol	Description
	Type BF
	Warning: See operation manual
<b>%SpO<sub>2</sub></b>	The pulse oxygen saturation
<b>♥ bpm</b>	Pulse rate (beats per minute)
	Low battery voltage
SN	Serial number
	Not for continuous monitoring

USA Corporate Headquarters:

**GF Health Products, Inc.**  
2935 Northeast Parkway  
Atlanta, Georgia 30360  
© 2009 GF Health Products, Inc.  
Made in China  
telephone: 800-347-5678, 770-447-1609  
fax: 800-726-0601, 678-291-3232

**GRAHAM-FIELD**  
[www.grahamfield.com](http://www.grahamfield.com)

Graham-Field and John Bunn are registered trademarks of GF Health Products, Inc.

Packaging, warranties, products, and specifications are subject to change without notice.

GF Health Products, Inc. is not responsible for typographical errors.

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Graham-Field y John Bunn son marcas registradas de GF Health Products, Inc.
GF Health Products, Inc. no es responsable por errores tipográficos.